


Medical Display Monitor RadiForce® RX320





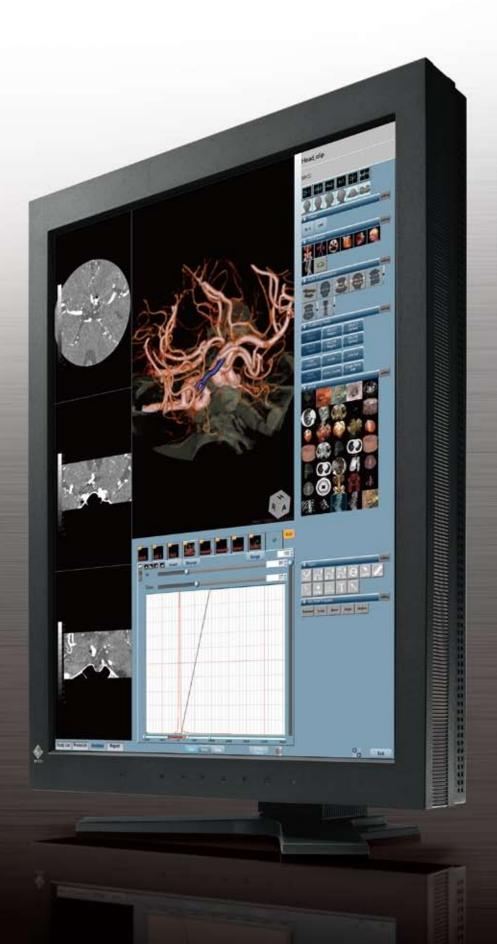

The high brightness 3MP color monitor is ideal for accurate display of chest X-ray, CR, and DR grayscale images and color images such as 3D color rendering and image fusion.

- High-quality IPS panel with 900 cd/m² brightness and 1000:1 contrast ratio which are on a par with the figures of monochrome monitors.
- 10-bit (RGB) color input compatibility which displays up-to 1.06 billion colors simultaneously providing highcrisp images for modality applications and 3D image fusion.
- Digital Uniformity Equalizer (DUE) function provides optimum brightness uniformity.
- CAL Switch function allows for various calibration modes to be selected from the front panel buttons for different modality images.
- **Dual link input support which reproduces a smooth moving picture without any dropping of frames.**
- Backlight Saver function reduces power consumption during prolonged usage.

Cabinet Color	Black	Power Management	DVI DMPM
Panel	TFT Color LCD Panel (IPS)	OSD Languages	English, French, German, Italian, Japanese, Simplified Chinese,
Active Display Size (H x V)	323.7 × 431.6 mm		Spanish, Swedish, Traditional Chinese
Viewable Image Size	Diagonal: 540 mm	Tilt / Swivel	40° Up / 35° Right, 35° Left
Native Resolution	1536 x 2048	Net Weight	With Stand: 10.4 kg, Without Stand: 7.4 kg
Pixel Pitch	0.21075 x 0.21075 mm	Hole Spacing	VESA standard 100 × 100 mm
Display Colors	10-bit colors: 1.06 billion (maximum) colors 8-bit colors: 16.77 million from a palette of 1.06 billion colors	Certifications and Standards	CE (Medical Device Directive), TÜV/GM (EN60601-1), cTÜVus (UL60601-1, CSA C22.2 No. 601-1), CB (IEC60601-1), VCCI-B, FCC-B, Canadian ICES-003-B, c-Tick, FDA 510(k)*, RoHS
Viewing Angles (H,V)	170° , 170°	Supplied Accessories AC power cord	AC power cord, user's manual, dual link signal cable (DVI-D
Brightness	900 cd/m ² (typical)		~ DVI-D), USB cable, EIZO LCD Utility Disk (RadiCS LE,
Recommended Brightness for Calibration	400 cd/m ²	Warranty	ScreenManager Pro for Medical), warranty card
Contrast Ratio	1000: 1 (typical)	Dimensions [Unit: mm]	Height Adjustment 82
On/Off Response Time	20 ms (typical)		
Scanning Frequency (H, V)	Digital: 31 - 127 kHz, 29 - 61 Hz (VGA Text: 69 - 71 Hz)		
	Frame synchronous mode: 59 - 61 Hz		
Dot Clock	215 MHz		
Input Signals	DVI Standard 1.0		
Input Terminal	DVI-D 24 pin		
USB Ports	1 upstream, 2 downstream		
USB Standard	Standard Rev. 2.0		
Power Requirements	AC 100 - 120 V, 200 - 240 V: 50 / 60 Hz		
Power Consumption	110 watts, Power Save Mode: Less than 3 watts		
Sensor	Backlight Sensor		

* Display of mammography images for diagnosis is not supported.



EIZO NANAO CORPORATION

153 Shimokashiwano, Hakusan, Ishikawa 924-8566 Japan Phone +81-76-277-6792 Fax +81-76-277-6793 radiforce.com All product names are trademarks or registered trademarks of their respective companies. RadiCS and RadiNET are trademarks, and EIZO, RadiForce, and ScreenManager are registered trademarks of Eizo Nanao Corporation. Specifications are subject to change without notice.